


AAJT-S APPLICATION & REMOVAL GUIDE

STEP-BY-STEP INSTRUCTIONS FOR SAFE AND EFFECTIVE USE




compression**works**

REMOVAL

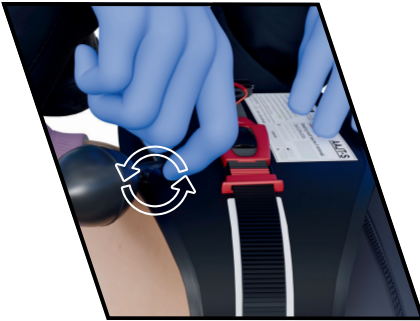

1 DEFLATE THE
BLADDER SLOWLY.

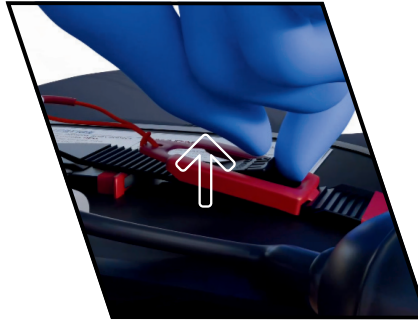

2 LIFT AND RELEASE
THE BLACK PULL TAB.



3 REMOVE THE
LADDER STRAP.

APPLICATION



1
CONNECT
RED TO RED.




2
POSITION
THE DEVICE



3
REMOVE SLACK
FROM BELT



4
RATCHET TO
TIGHTEN

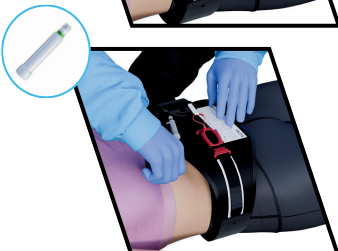

5
INFLATE BLADDER
UNTIL **GREEN**

DO NOT CUT



PREMATURE REMOVAL
MAY RESULT IN
RESUMED BLEEDING



SCAN THE QR CODE
TO WATCH A
HOW-TO-USE VIDEO.

WWW.COMPRESSIONWORKS.COM